




GeneProof®

GeneProof Mycoplasma pneumoniae PCR Kit

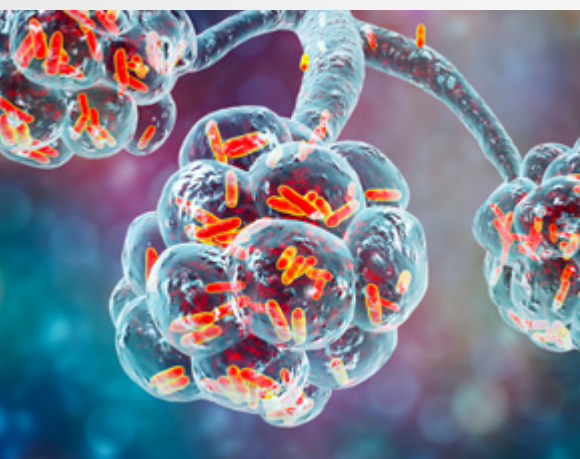

HIGH SENSITIVITY AND SPECIFICITY

- Secured by targeting the multicopy sequence of the *M181* gene

CONTAMINATION PREVENTION

- Master Mix contains Uracil-DNA glycosylase (UNG) and dUTPs eliminating carryover contamination

EASY-TO-USE CONCEPT

- Single tube Ready-to-Use Master Mix contains all components for PCR amplification
- No additional PCR reagents pipetting necessary

SIMPLE LABORATORY WORKFLOW

- Easily combinable with other GeneProof PCR kits in one workflow

ORDER INFORMATION

REF	PACKAGE
MP/ISEX/025	25 reactions
MP/ISEX/100	100 reactions



COMPATIBLE WITH A WIDE
RANGE OF REAL-TIME PCR
DEVICES



CERTIFIED
DIAGNOSTIC TEST



GeneProof *Mycoplasma pneumoniae* PCR Kit

- + GeneProof Adenovirus PCR Kit
- + GeneProof Aspergillus PCR Kit
- + GeneProof Bordetella pertussis/
parapertussis PCR Kit
- + GeneProof Enterovirus PCR Kit

- + GeneProof Chlamydia
pneumoniae PCR Kit
- + GeneProof Flu Multiplex PCR Kit
- + GeneProof Legionella
pneumophila PCR Kit

- + GeneProof *Mycoplasma pneumoniae* PCR Kit
- + GeneProof Mycobacterium
tuberculosis PCR Kit

GENEPROOF COVID-19 SOLUTION

- + GeneProof SARS-CoV-2
PCR Kit

- + GeneProof SARS-CoV-2
Screening PCR Kit

- + GeneProof SARS-CoV-2
Advanced PCR Kit

INDICATION	<i>in vitro</i> diagnostic medical device
REGULATORY STATUS	CE IVD
INTENDED USER	For professional use in laboratories with trained staff
TECHNOLOGY	Real-time PCR
TYPE OF ANALYSIS	Qualitative
TARGET SEQUENCE	The <i>M181</i> gene encoding the CARDS toxin
ANALYTICAL SPECIFICITY	<i>M. pneumoniae</i> , 100%
ANALYTICAL SENSITIVITY (LOD)	0.46 cp/μl
EXTRACTION/INHIBITION CONTROL	PCR inhibition and DNA extraction efficiency control
VALIDATED SPECIMEN	BAL, sputum, swab
STORAGE	-20 ± 5 °C
VALIDATED EXTRACTION METHODS	croBEE 201A Nucleic Acid Extraction System GeneProof PathogenFree DNA Isolation Kit
INSTRUMENTS	croBEE Real-Time PCR System Applied Biosystems 7300 / 7500 Real-Time PCR System AriaMx Real-Time PCR System AMPLilab Real-Time PCR System BioQuant-96 Real-Time PCR System CFX Connect™ / CFX96™ / Dx Real-Time PCR Detection System Gentier 96E/96R Real-Time PCR System LightCycler® 2.0 / 480 LineGene 9600 / 9600 Plus Mic qPCR Cyclor QuantStudio™ 3 / 5 Real-Time PCR System Rotor-Gene 3000 / 6000 / Q SLAN® Real-Time PCR System StepOne™ / StepOnePlus™ Real-Time PCR System
REQUIRED DETECTION CHANNELS	FAM, HEX/JOE/VIC
EXTERNAL QUALITY ASSESSMENT	Regularly tested in QCMD and Instand e.V. External Quality Assessment Panels – results at www.geneproof.com